

HidroDent® Oral Irrigator



ORAL IRRIGATOR

Model NO.: HD-3000



IMPORTANT SAFEGUARDS



When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:
READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not place charger in or drop charger into water or other liquid.
- Check the charger cord for damage before the first use and during the life of the product.


IMPORTANT SAFEGUARDS WARNING: To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- If product is opened/disassembled for ANY reason, warranty is VOID.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Do not direct water under the tongue, into the ear, eyes, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet.
- Only use tips and accessories that are recommended by HydroDent
- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Close supervision is necessary when this product is used by, on, or near children or invalids.


SAVE THESE INSTRUCTIONS.



PRODUCT DESCRIPTION

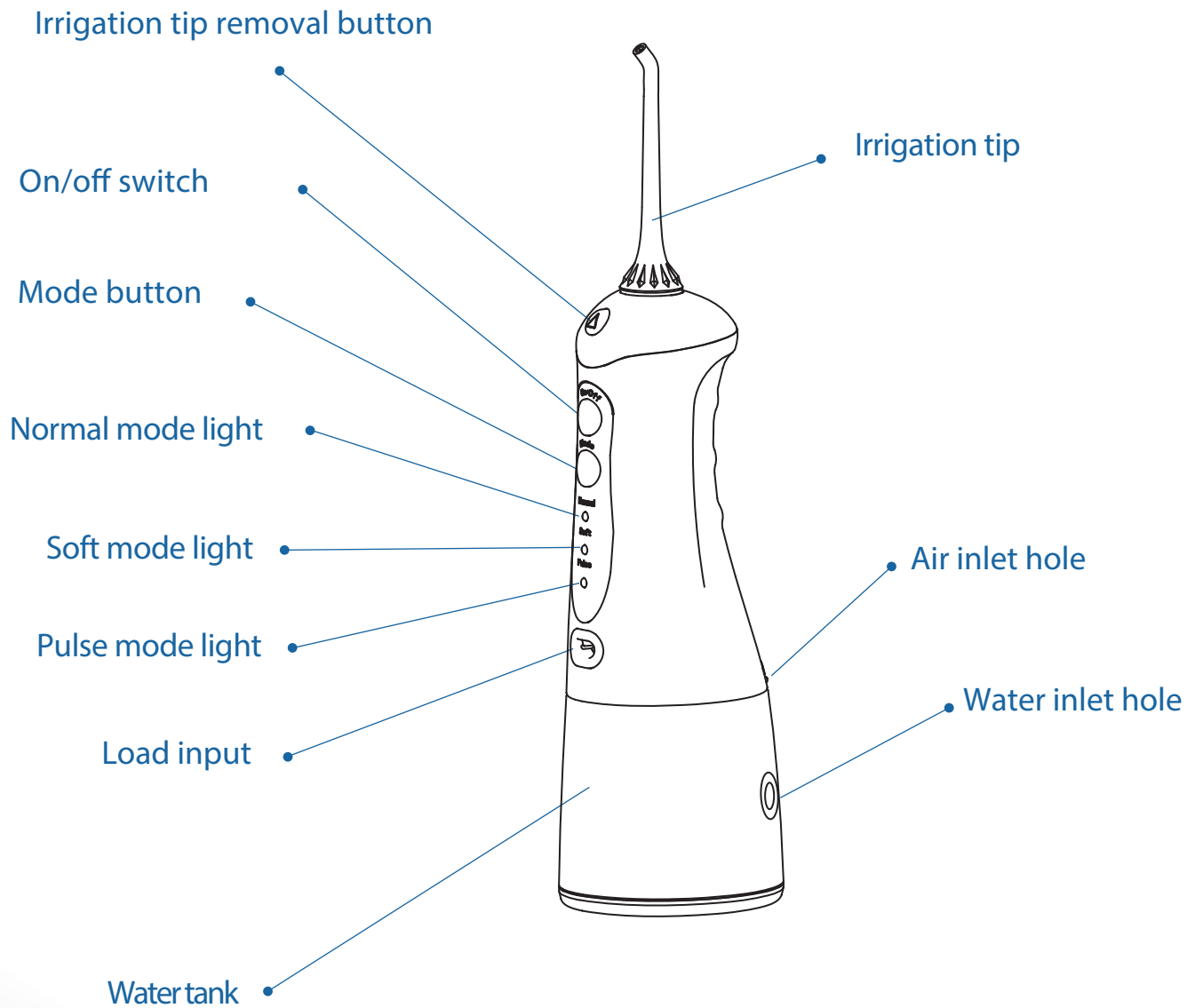


PRODUCT ACCESSORIES

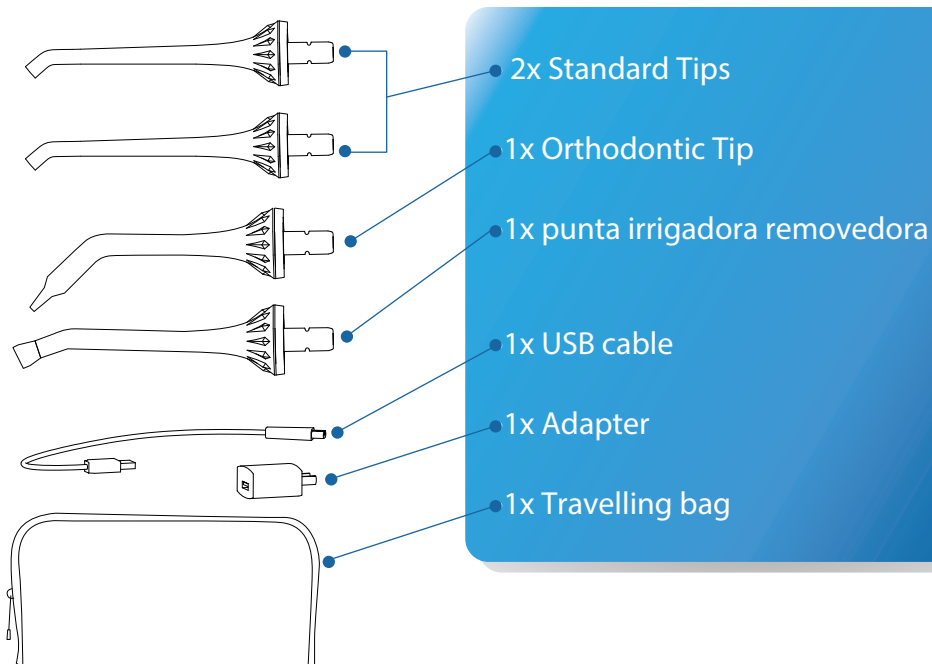

Accessories can be purchased in www.hidrodent.com

OPERATING PROCEDURES

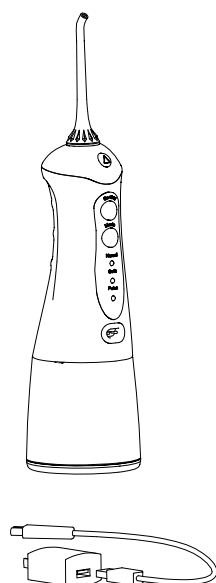


CHARGING UNIT PRIOR TO FIRST USE


Plug the power cord into the wall and into the HidroDent® for approximately 4 hours. The lights will then turn off when the unit is fully charged. Although it may turn off before 4 hs, we recommend you keep it plugged for at least 4 hours. A full charge will last for approximately 15 to 20 days of typical use. When the batteries need of charging, the HidroDent® will have less pressure, and will stop working. A complete recharge will take 4 hours.




INSTALLATION AND DISASSEMBLY OF TIP

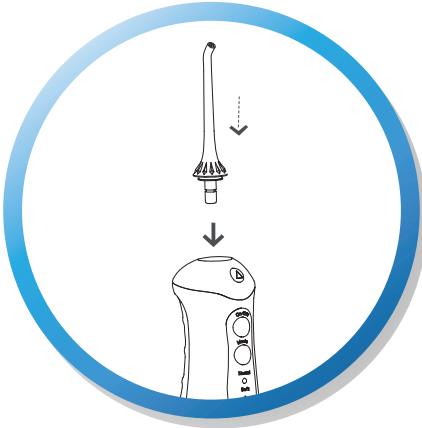

Insert the jet tip into the spray head hole until you hear a “click”

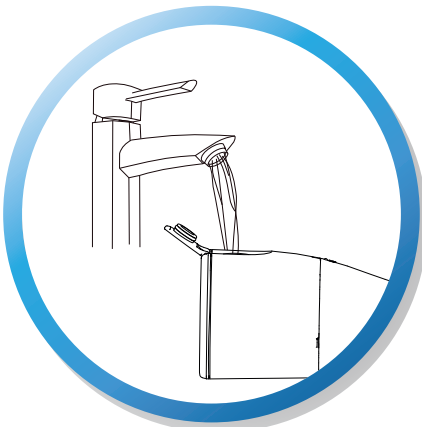


Open the water tank and pour warm water into the water tank.
WARNING: Please ensure that when pouring water into the water tank, the device is switched off.

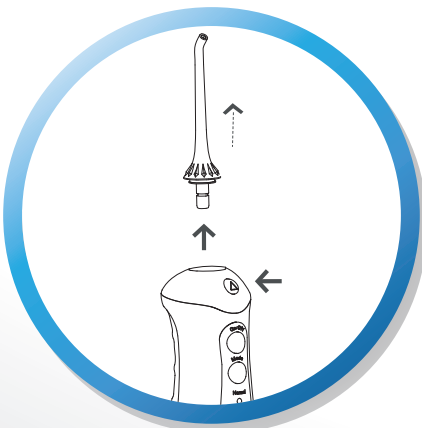

When removing or replacing the irrigator tip, please press the release button on the handle of the irrigator tip before removing it.




- This product has a memory function, and after switching it back on the water pressure mode will be the same as when the device was switched off.
- The normal, soft and pulse modes may be switched whether the device is on or off.
- The modes may be changed when the device is switched off, and you can select the mode you like before switching the device back on.
- Irrespective of the mode, after operating for 120 seconds the product will automatically switch off.
- During recharge, the three display lights will flash. When all three display lights are constantly on, this indicates that the power recharged.
- During usage, if three display light are flashing, this indicates that the product's power is insufficient, and recharging is required.
- During recharge this product cannot be used.
- Rotate the jet tip dial in order to adjust the direction of water flow

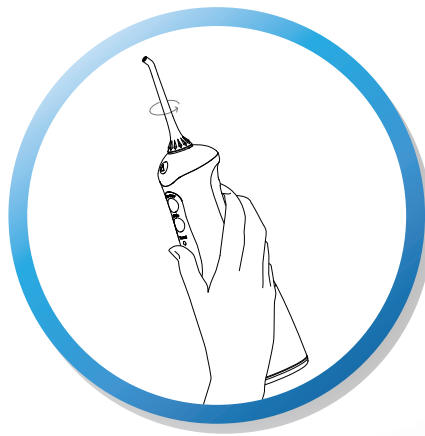





RECOMMENDED TECHNIQUES

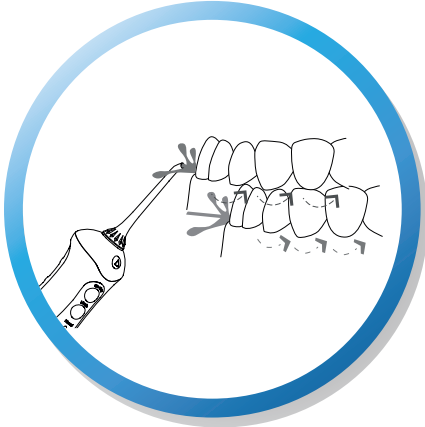

Holding the device at downward slant, position the jet tip into the mouth, and direct the jet tip at teeth before switching on the device.

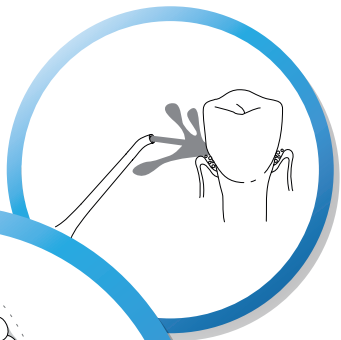


When the water flow and the gums are at a 90 degree angle, slightly close your mouth in order to avoid the water spraying out, and achieving the best results. Please first clench the teeth (rear teeth) and begin to apply to the incisors. Slide the water flow along the gums, as well as pause briefly between teeth, until the area around the gums and the spaces between the teeth are fully cleaned.

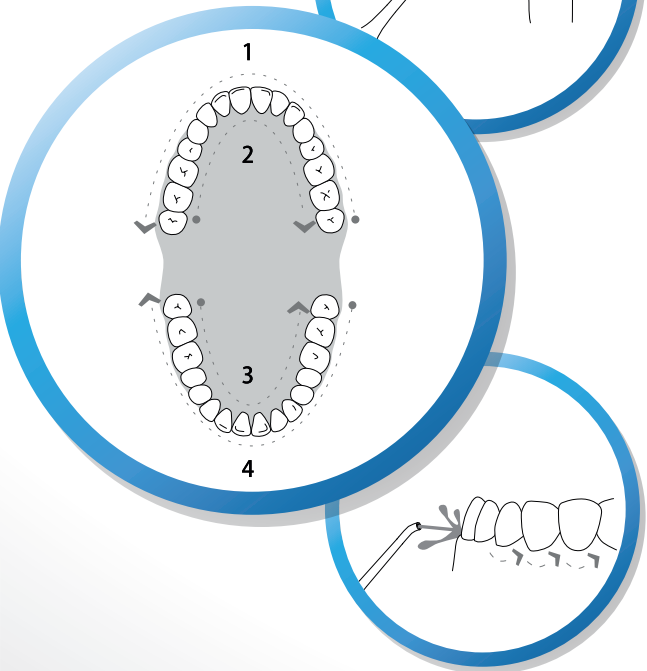

STEPS TO FOLLOW


1

Turn off the device and empty the water tank

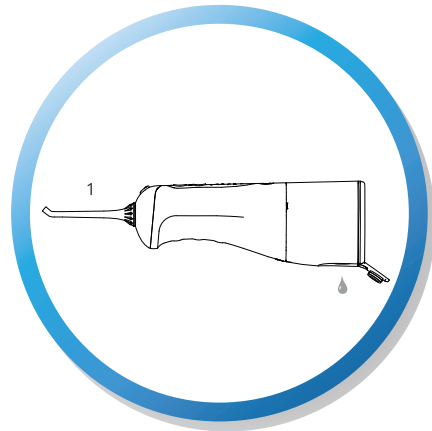

2

Using both hands as indicate in the Diagram – Use your right hand to hold the water tank. At the same time rotate downwards in an anti-clockwise direction, in order to remove the water tank.

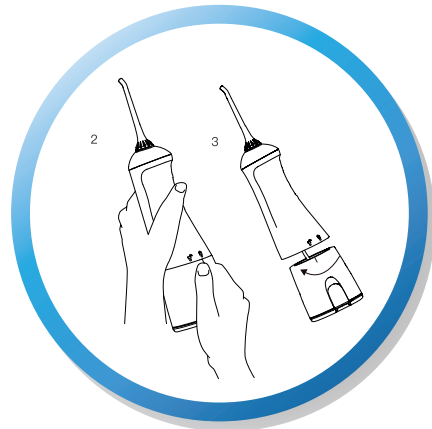



USE AND MAINTENANCE OF THE BATTERY

BATTERY MAINTENANCE

- Prior to initial usage the device must be fully charged, for approximately 4 hours.
- If to be stored for a long-period of time, please fully charge the device in advance.
- In order to maximize the usage life of the battery, please promptly recharge when power is insufficient.

BATTERY USE

The device includes a lithium battery which cannot be changed. Towards the end of the battery life, please dispose in a local recycling centre



Problems	Reasons	Method to handling
Leakage of water between jet tip and screw part	The seal ring on the fixed socket for the jet tip is damaged, and the jet tip cannot be installed properly	Contact the retailer for repair, and install the jet tip properly
Water tank leaks water	Seal ring is loose	Re-install correctly
Inadequate pressure	Low voltage usage	Use after recharging
Device cannot be switched on	Low voltage usage	Use after recharging
Product is not discharging water	Filtration net is clogged	Clean the filtration net
Water tank sides are leaking water	Air entry hole is producing water	Standard phenomenon – high air pressure principle



If the above methods cannot resolve these problems, please consult with the retailer or proceed to the designated after-sales repair site for repairs.

Gentle reminder:

1. At the time of initial usage, please select “Soft mode”, in order to avoid the gums bleeding due lack of adaption to usage. After consistent usage for weeks there should be improvements. If bleeding is excessive or bleeding continues after one week, Chronic gum inflammation could be an issue, and it is recommended that professional dental treatment be sought.
2. It's not recommended that the product be used by children, or people with dental implants, oral ulcers or open oral sores.

SERVICES POST-SALE

Thank You for Purchasing The HidroDent Portable Oral Irrigator – Please Take Note of the Following Issues:



1 Please retain the effective sales invoice issued by retailer as a maintenance document.

2 One year product malfunction warranty for non-human malfunctions is provided for this product accordance with the relevant state provisions. If the product was purchased via HidroDent.com, the documentation is not required, and only the order form number needs to be provided, in the order form number needs to be provided, in order for us to provide corresponding after-sales free maintenance (based on the invoice date for purchase of the device).

3 Where the following circumstances arise during the maintenance period, cost and services fee shall be charged based on circumstances;

- The customer causes damage due to inappropriate usage or maintenance;
- Repairers who are not covered by third-party repair cause damage during repairs;
- Lack of a relevant sales document or more than two years have elapsed since the date of manufacture;
- Damage caused during usage for commercial purposes;
- Damage caused by force majeure events (such as storms, earthquake or fire).



ONE YEAR LIMITED WARRANTY

Maintenance service

HidroDent® irrigators have no consumer service able electrical components and do not require routine service maintenance. For all parts and accessories, visit www.hidrodent.com

You can also call our customer service number at +507 3086967 (Panamá), or you can write to us at info@hidrodent.com

Still have questions? Visit www.hidrodent.com. HidroDent® guarantees the Original consumer of this new product that is free from defects in material and workmanship for 1 year from the date of purchase. The consumer must submit the original purchase receipt as proof of the purchase date and, if requested, send us the entire product to support a warranty claim. HidroDent® will replace at its discretion any part of the product, provided that the product has not been abused, misused, altered, or damaged after purchase, was used in accordance with the instructions, and used only with accessories or consumable parts approved by HidroDent®.

This limited warranty excludes accessories or consumable parts such as tips, etc. Use with the wrong voltage will damage the product and void the warranty. This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.