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HydroDent

Hydro flosser system
Model HD-3000 series

In order to facilitate the warranty process, and receive incredible benefits, please register your product here and forget about the receipt: hidrodent.com/registro


CONGRATULATIONS!

You have just acquired the perfect complement to brushing, and you have started the path towards a healthier smile.


At HydroDent®, your dental health is the most important. We want our clients to be 100% satisfied with the product and the results obtained. If you have any issue, PLEASE CONTACT US BEFORE RETURNING THE PRODUCT to, info@hidrodent.com or contact us through our website: www.hidrodent.com

Quick Start Guide


- 1** Charge the unit for 4hs, connecting it to the charger.

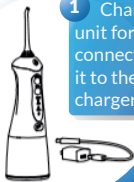- 2** Insert the tip of your choice until it clicks. To eject it, press the button. Do not use without a tip, and do not eject the tip while running

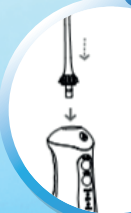- 3** Open the lid of the tank and fill it with warm water. You can add mouthwash if you want.

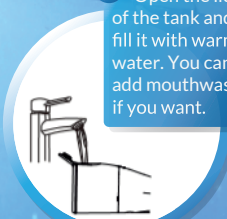
- 4** Press the "mode" button and adjust to the desired pressure. Start with soft pressure until you get used to it, and then increase as desired.

- 5** Insert the tip into your mouth, close your lips to avoid splashing, and then press the ON button. Try to keep the HydroDent in a 90° position to ensure continuous flow.

- 6** Start by aiming at the gum, and run the entire line through each tooth, pausing briefly between each tooth.

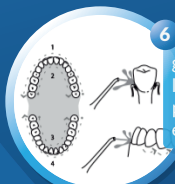- 7** Once irrigation is complete, empty the tank completely, and leave it open to avoid humidity and condensation inside.

- 8** To rinse and clean the tank and keep the equipment in perfect condition, dry the HydroDent on the outside, and holding it steady with one hand, turn the tank counterclockwise and remove it gently.

TO AVOID SPLASHING, KEEP IT INSIDE YOUR MOUTH WITH ENCLOSED LIPS WHILE THE UNIT IS ON.

IT MAY BLEED IN THE GUMS DURING THE FIRST USES, DON'T PANIC, IT WILL DISAPPEAR AFTER FEW USES. IF BLEEDING PERSISTS FOR MORE THAN A WEEK, CONSULT YOUR DENTIST.